

Extremely user-friendly

- Simple control panel
- Big cutting table
- Cleaning gun with strong water flow
- A powerful LED lamp lights up the cutting area



Labotom-5 is a high-quality table-top manual cut-off machine for 250 mm / 10" dia. cut-off wheels. It is specially designed to be used both in labs and production areas.

With its simple control panel, the ergonomic handle as well as safe lock and emergency stop, the Labotom-5 is not only extremely user-friendly - it also provides optimal protection for its users.

The large cutting table offers extra space to position your workpiece and after working on the Labotom-5, the efficient flushing of the chamber makes cleaning up fast and simple.

Safety lock and emergency stop

Labotom-5 is equipped with both a safety lock and an emergency stop for optimal protection of the users. The motor is braked electrically which ensures that it is stopped before the door safety lock is released.




Cleaning gun with strong water flow

For easy maintenance Labotom-5 is equipped with a cleaning gun. The cleaning gun is easily accessible and stored in a holder on the side of the machine.

Dedicated cut-off wheels


A complete range of dedicated Struers' cut-off wheels is available for Labotom-5. Struers cut-off wheels are of the highest quality and maintain homogeneous cutting characteristics throughout their life, due to built-in wear compensation and thereby provides better control of the cut.

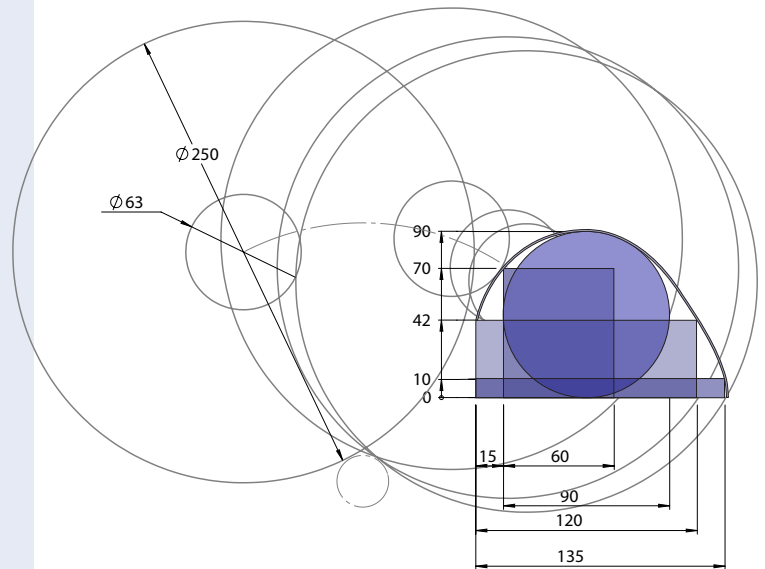

Cutting capacity

The projected cutting capacity for a new cut-off wheel is up to 90 mm / 3.5" dia. or 120 x 42 mm / 4.7" x 1.7". Please refer to the illustration.

To cut workpieces exceeding the width of the cutting chamber in the left side, a plate in the cover may be removed and replaced by a Safety Extension Box.

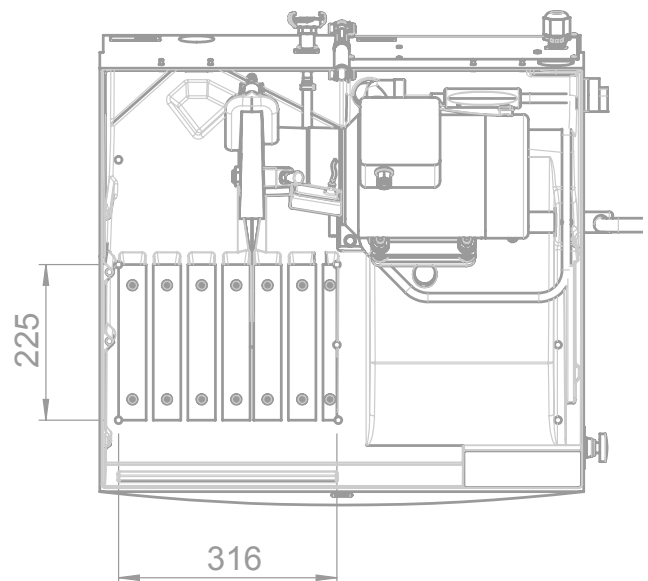

Big cutting table


Labotom-5 has a big cutting table (316 mm / 12.4" wide and 225 mm / 8.9" dept) for easy clamping of the workpiece. The cutting table has five 10 mm T-slots; three on the left hand side and two on the right hand side. The big cutting table ensures that almost all shapes of workpieces can be clamped and cut safely.

Big slope in the cutting chamber

The sloping base of the cutting chamber prevents cooling water from accumulating in the chamber. A prism knob is placed in the drainage hole so that cut-off specimens will not enter the drain.

Accessories

Cat. no.:

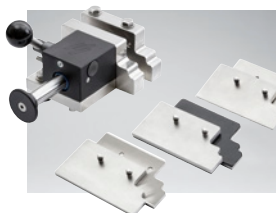

Quick Clamping Device

For securing the workpiece. To be mounted on the left side of the cut-off wheel.
Complete with backstop. Max clamping width: 50 mm.

06316912



Vertical Clamping Device

For clamping irregularly shaped workpieces on the cutting table. Clamping height is adjustable up to 58 mm. Complete with operating key and one flat clamping shoe. Used for 10 mm T-slot.

05116905

Swivel shoes for Vertical Clamping Device

Set of 4 multi-shaped swivel shoes for Vertical Clamping Device (05046904, 05116905 or 05036915).

05116911

Arm Extension

For extending arm of Vertical Clamping Device (05036915, 05046904 and 05116905)

05046906



Vertical Quick Clamping Device

For securing the workpiece. To be mounted on the left side of the cut-off wheel.
Max clamping width: 50 mm

06316922



Single Prism for Vertical Quick Clamping Device

Jaw block for optimum clamping of cylindrical and oval workpieces.
For use with Vertical Quick Clamping Device for 10 mm T-slots (05876925).

05876926



Spring Clamp

For retaining the workpiece. Complete with backstop. To be mounted on Discotom-60/-65, Discotom-6 and Labotom-5 at the right side of the cut-off wheel. For 10 mm T-slots.

05136903



Fixed Stand

Stand for specimen holders with dovetail clamp. 40 mm manual positioning. For 10 mm T-slot.

05876927



Specimen Holder CATFF 04946904

Teardrop type with opening.
For workpieces 15-40 mm dia.



Specimen Holder CATTY

With 7 screws for irregular workpieces.
Max 40 mm width.

04946905



Specimen Holder CATSI 04946903

For general use.
Max opening: 60 mm.



Specimen Holder CATRE

For cylindrical workpieces.
Max opening 40 mm dia.

04946907



Safety Extension Box

For extending the cutting chamber when cutting long workpieces.
To be mounted on the left side of Labotom-5.
For workpieces with a maximum size of
H x W: 90 mm / 3.5" dia. or 120 x 42 mm / 4.7" x 1.7"
Dimensions, W x D x H: 652 x 120 x 95 mm / 6.5" x 4.7" x 3.7"

06046903

For more information regarding clamping tools and cut-off wheels please see special brochures.

Technical data

Cut-off wheel

Size	250 mm / 10" dia.
Arbor size	32 mm / 1.26" dia.
Rotational speed (running idle)	2845 rpm (50 Hz) / 3450 rpm (60 Hz)
Cutting capacity *	90 mm dia. / 3.5" dia. or 42 x 120 mm or 1.7" x 4.7"

Cutting table

Width	316 mm / 12.4"
Depth	225 mm / 8.9"
T-slots	10 mm

Cutting motor

Voltages	Nominal load S3 (60%)	Intermittent S3 (15 %)	Fuse
3 x 200 V / 50 Hz	2.2 kW / 3.0 HP	3.2 kW / 4.3 HP	40 A
3 x 200 - 210 V / 60 Hz	2.6 kW / 3.5 HP	3.2 kW / 4.3 HP	30 A
3 x 220 - 230 V / 50 Hz	2.2 kW / 3.0 HP	3.2 kW / 4.3 HP	40 A
3 x 220 - 240 V / 60 Hz	2.6 kW / 3.5 HP	3.2 kW / 4.3 HP	30 A
3 x 380 - 415 V / 50 Hz	2.2 kW / 3.0 HP	3.2 kW / 4.3 HP	40 A
3 x 380 - 415 V / 60 Hz	2.6 kW / 3.5 HP	3.2 kW / 4.3 HP	30 A
3 x 460 - 480 V / 60 Hz	2.6 kW / 3.5 HP	3.2 kW / 4.3 HP	30 A

Dimensions

Width	700 mm / 27.6"
Width with handle and liquid gun	960 mm / 37.8"
Depth	760 mm / 29.9"
Depth with plug	808 mm / 31.8"
Height, machine	440 mm / 17.3"
Height, handle	520 mm / 20.5"
Height, open cover	942 mm / 37.1"
Weight	75 kg / 165 lbs

Surrounding temperature	5 - 40 °C / 41-104 °F
--------------------------------	-----------------------

Humidity	35 - 85 % RH, non condensing
-----------------	------------------------------

* The actual cutting capacity depends on the specific workpiece dimensions, cut-off wheel and clamping technique.

Specifications

Labotom-5

Manual cut-off machine. Motor power 2.2 / 2.6 kW (3.0 / 3.5 HP). For 250 mm (10") cut-off wheels.

With cutting table with 10 mm T-slots. Recirculation cooling unit (05766xxx) and clamping tools are ordered separately. Optional table unit (06266101) can be used. A selection guide for cut-off wheels is included in the digital consumables catalogue on the Struers website.

Cat. no.:	Voltage
06046229	3 x 200 V / 50 Hz
06046230	3 x 200-210 V / 60 Hz
06046235	3 x 220-230 V / 50 Hz
06046236	3 x 220-240 V / 60 Hz
06046246	3 x 380-415 V / 50 Hz
06046247	3 x 380-415 V / 60 Hz
06046254	3 x 460-480 V / 60 Hz

Struers' products are subject to constant product development.
Therefore, we reserve the right to introduce changes in our products without notice.

AUSTRALIA & NEW ZEALAND
Struers Australia
27 Mayneview Street
Milton QLD 4064
Australia
Phone +61 7 3512 9600
Fax +61 7 3369 8200
info.au@struers.dk

BELGIUM (Wallonie)
Struers S.A.S.
370, rue du Marché Rollay
F- 94507 Champigny
sur Marne Cedex
Téléphone +33 1 5509 1430
Télécopie +33 1 5509 1449
struers@struers.fr

BELGIUM (Flanders)
Struers GmbH Nederland
Zomerdijk 34 A
3143 CT Maassluis
Telefoon +31 (10) 599 7209
Fax +31 (10) 5997201
netherlands@struers.de

CANADA
Struers Ltd.
7275 West Credit Avenue
Mississauga, Ontario L5N 5M9
Phone +1 905-814-8855
Fax +1 905-814-1440
info@struers.com

CHINA
Struers Ltd.
No. 1696 Zhang Heng Road
Zhang Jiang Hi-Tech Park
Shanghai 201203, P.R. China
Phone +86 (21) 6035 3900
Fax +86 (21) 6035 3999
struers@struers.cn

CZECH REPUBLIC & SLOVAKIA
Struers GmbH Organizační složka
vědeckotechnický park
Přílepká 1920,
CZ-252 63 Roztoky u Prahy
Phone +420 233 312 625
Fax +420 233 312 640
czechrepublic@struers.de
slovakia@struers.de

GERMANY
Struers GmbH
Carl-Friedrich-Benz-Straße 5
D- 47877 Willich
Telefon +49 (0) 2154 486-0
Fax +49 (0) 2154 486-222
verkauf@struers.de

FRANCE
Struers S.A.S.
370, rue du Marché Rollay
F-94507 Champigny
sur Marne Cedex
Téléphone +33 1 5509 1430
Télécopie +33 1 5509 1449
struers@struers.fr

HUNGARY
Struers GmbH
Magyarország Fióktelepe
2040 Budaörs
Szabadság utca 117
Phone +36 2380 6090
Fax +36 2380 6091
Email: hungary@struers.de

IRELAND
Struers Ltd.
Unit 11 Evolution@ AMP
Whittle Way, Catcliffe
Rotherham S60 5BL
Tel. +44 0845 604 6664
Fax +44 0845 604 6651
info@struers.co.uk

ITALY
Struers Italia
Via Monte Grappa 80/4
20020 Arese (MI)
Tel. +39-02/38236281
Fax +39-02/38236274
struers.it@struers.it

JAPAN
Marumoto Struers K.K.
Takanawa Muse Bldg. 1F
3-14-13 Higashi-Gotanda,
Shinagawa
Tokyo
141-0022 Japan
Phone +81 3 5488 6207
Fax +81 3 5488 6237
struers@struers.co.jp

Struers ApS

Pederstrupvej 84
DK-2750 Ballerup, Denmark
Phone +45 44 600 800
Fax +45 44 600 801
struers@struers.dk
www.struers.com

NETHERLANDS
Struers GmbH Nederland
Zomerdijk 34 A
3143 CT Maassluis
Telefoon +31 (10) 599 7209
Fax +31 (10) 5997201
netherlands@struers.de

NORWAY
Struers ApS, Norge
Sjøslogenveien 44C
1407 Vinterbro
Telefon +47 970 94 285
info@struers.no

AUSTRIA
Struers GmbH
Zweigniederlassung Österreich
Betriebsgebiet Puch Nord 8
5412 Puch
Telefon +43 6245 70567
Fax +43 6245 70567-78
austria@struers.de

POLAND
Struers Sp. z o.o.
Oddział w Polsce
ul. Jasnogórska 44
31-358 Kraków
Phone +48 12 661 20 60
Fax +48 12 626 01 46
poland@struers.de

ROMANIA
Struers GmbH, Sucursala Bucuresti
Str. Preciziei nr. 6R
062203 sector 6, Bucuresti
Phone +40 (31) 101 9548
Fax +40 (31) 101 9549
romania@struers.de

SWITZERLAND
Struers GmbH
Zweigniederlassung Schweiz
Weissenbrunnstraße 41
CH-8903 Birmensdorf
Telefon +41 44 777 63 07
Fax +41 44 777 63 09
switzerland@struers.de

SINGAPORE
Struers Singapore
627A Aljunied Road,
#07-08 BizTech Centre
Singapore 389842
Phone +65 6299 2268
Fax +65 6299 2661
struers.sg@struers.dk

SPAIN
Struers España
Camino Cerro de los Gamos 1
Building 1 - Pozuelo de Alarcón
CP 28224 Madrid
Teléfono +34 917 901 204
Fax +34 917 901 112
struers.es@struers.es

FINLAND
Struers ApS, Suomi
Hietalahdenranta 13
00180 Helsinki
Puhelin +358 (0)207 919 430
Faksi +358 (0)207 919 431
finland@struers.fi

SWEDEN
Struers Sverige
Box 20038
161 02 Bromma
Telefon +46 (0)8 447 53 90
Telefax +46 (0)8 447 53 99
info@struers.se

UNITED KINGDOM
Struers Ltd.
Unit 11 Evolution @ AMP
Whittle Way, Catcliffe
Rotherham S60 5BL
Tel. +44 0845 604 6664
Fax +44 0845 604 6651
info@struers.co.uk

USA
Struers Inc.
24766 Detroit Road
Westlake, OH 44145-1598
Phone +1 440 871 0071
Fax +1 440 871 8188
info@struers.com