
HDT AND VICAT TESTING SYSTEMS

Introducing the HV Series Testing Systems




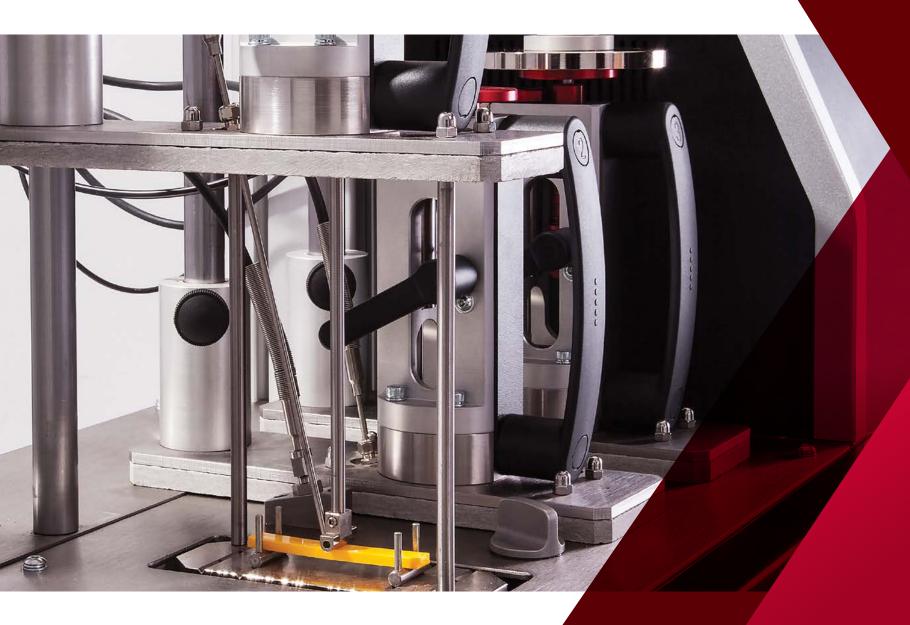

Instron brings over 50 years of experience in plastics and composites testing solutions to thousands of customers globally. Supported by a world class service team, we deliver the best solutions to laboratories that are looking to grow their business and efficiently verify the quality of their products.

INTRODUCING THE HV SERIES

The design of the HV Systems and BLUEHILL® HV Software is based on the full understanding of customers needs for both an enhanced user-experience and outstanding productivity, making materials testing SIMPLER, SMARTER, SAFER.

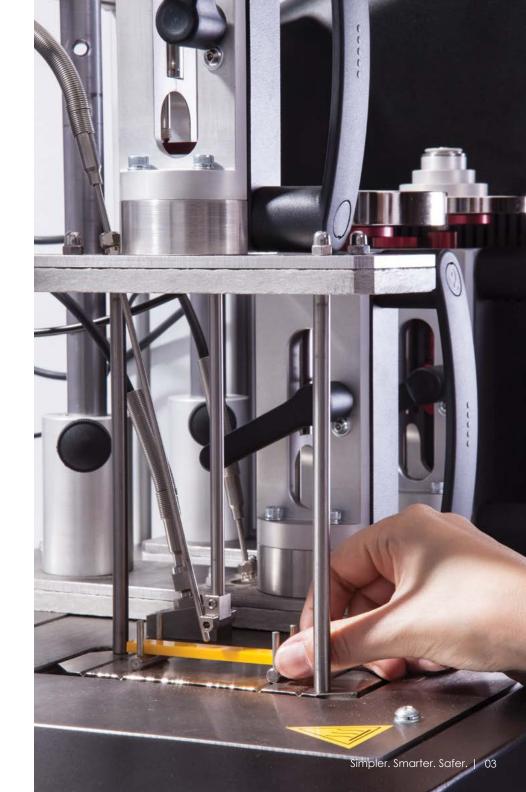

SIMPLER Fast and Effortless Workflow





The intuitive design of weights eliminates the risk of selection mistakes. The weights are stackable and enable the user to handle them with confidence. A wide range of stress values can be achieved without the need for any additional tools.



Comfortable Handling of the Testing Stations

The smooth lifting mechanism of the stations, combined with the ergonomic position of the handle, allows the user of the HV3S to raise and lower the stations safely and easily.





Touch Screen Interface

The user-friendly touch screen has a simple and clean look and allows for operation without a computer. The workflow is intuitive and offers a seamless interaction with the system for a quick test setup and data export. The pre-formatted templates enalble test run in a few simple touches.

Effortless Workflow

The optional Bluehill HV Software enhances the ease of use. The user can run tests, edit methods, analyze results or configure the system with a few clicks. The simplicity of the software allows the user to run tests with minimal effort increasing the laboratory efficiency and throughput.

SMARTER

Integrated Innovations for Efficient Tests



Oil degradation control



The system constantly monitors the oil, alerts when the bath level is out of range and reports when oil degradation that could adversely affect test results has occurred. This feature helps reduce the risk of running invalid tests and increases productivity. It's a perfect solution to maintain recommended oil bath conditions and extend the life of the machine.

Advanced Electronic Technology



The advanced electronics on the HV Series automatically zero the position of the LVDT measuring sensors before starting a test, reducing test time and human errors. Correction factors are calculated and applied to compensate for the thermal expansion at each station ensuring accurate results.

Efficient Cooling



The HV Series is designed for superior cooling efficiency, rapidly returning the oil to the starting temperature after a test run. This allows the user to begin the next test with minimal wait time and reduces water consumption costs.

Small Footprint




As a result of increasingly limited space in laboratories, equipment is often squeezed into tight spaces, creating crowded and less functional laboratories. The HV Series compact design helps maximize valuable laboratory space resulting in enhanced testing efficiency.

Superior Results and Repeatability



Factory trained engineers provide comprehensive verifications, checking the dimensions of critical components and temperature behavior, to assure results are accurate and repeatable. The verification certificates provide full traceability, providing confidence in the laboratory's results.

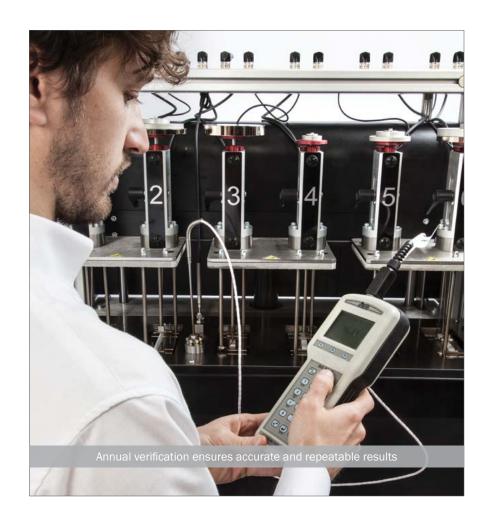

SAFER

Safety without Sacrificing Throughput







The full shield on the HV6X model restricts access during test execution, protecting the user from moving parts. It is the ideal solution for open access workplaces, where many testing instruments are simultaneously in use and a risk of test interruption exists.

Maximum Temperature Control



The HV Series has a built-in temperature control system which limits the maximum temperature allowed for the oil bath. This helps conform to the specifications of each oil type, and prevents improper temperature settings that can generate dangerous fumes and inaccurate test results.

Secured Test Access



Bluehill HV Software secures the test methods, preventing inadvertent changes, to ensure repeatability and standards compliance, over multiple shifts and numerous operators. The software suite allows for creation of specific user profiles based on the operator's level of expertise. Bluehill HV is tailored to simplify the testing process and ensure data integrity.



Effective Training Program

Instron engineers offer comprehensive on-site training courses to minimize operators' exposure to risk, improve efficiency, optimize the system's capabilities and to support the laboratory's health and safety policies to enhance the ownership experience.



SIMPLER. SMARTER. SAFER.

The HV Series Testing Systems

HV3S

The HV3S is a great choice for laboratories that perform HDT or Vicat tests at a low volume. The compact size and footprint make it an ideal system for busy laboratories with a lot of test equipment.

Manual lifting station and weights application | up to 3 stations

The HV SERIES offers two different bath dimensions: small to host up to 2 testing stations and large for up to 6 stations.



HV6M

Motorized lifting station and manual weights application | up to 6 stations

The HV6M is suitable for laboratories where the needing of HDT or Vicat testing is part of the daily activity. The motorized lifting mechanism guarantees less operator effort.



HV6X

Motorized lifting station and weights application | up to 6 stations | Interlocked cover shield

The HV6X is suitable for laboratories where there is a continuous backlog of specimens to be tested and reduced time to test is essential. In addition, it fits perfectly in those work environments where safety procedures are extremely strict and system automation is necessary.





RUN A TEST IN A FEW STEPS

Laboratories receive various materials that need to be tested on a daily basis. Bluehill HV has touch tiles with easy access to test run, method or analysis. The Software is developed around a simple and customizable workflow that systematically guides the user through a test, minimizing set up time. Bluehill HV improves productivity with its analysis module and helps perform testing in an intuitive touch-friendly way.





Test Method Compliance

Create ISO and ASTM test templates with smart default values so the user is ready to test with minimal setup time, increasing laboratory efficiency. Bluehill HV Software provides the ability to test samples to different methods simultaneously, helping to optimize throughput.



Analysis of Test Results

Test data, graphs and inputs are shown in a single screen allowing the user to analyze results at a glance. Set acceptance criteria to automatically reject out of tolerance samples. Create report templates with minimal effort. Export and print test reports or send test data to LIMS or TrendTracker.

SPECIFICATION TABLE


Simpler. Smarter. Safer







		HV3S	HV6M	HV6X
Number of testing stations	-	Up to 3 (includes 2)	Up to 6 (includes 2)	Up to 6 (includes 2)
Movement of the stations	-	Manual	Machine controlled	Machine controlled
Application of test weights	-	Manual	Manual	Machine controlled
Heat transfer medium	-	Silicone oil	Silicone oil	Silicone oil
Oil bath capacity	litres USA gal	9 2.4	16 4.2	16 4.2
emperature range	°C °F	20 - 300 68 - 572	20 - 300 68 - 572	20 - 300 68 - 572
emperature resolution	°C °F	0.1 0.18	0.1 0.18	0.1 0.18
VDT resolution	mm in	0.001 0.00004	0.001 0.00004	0.001 0.00004
imensions (W × D × H)	mm in	938 x 480 x 713 36.9 x 18.9 x 28.1	1291 x 480 x 780 50.8 x 18.9 x 30.7	1291 x 550 x 900 50.8 x 21.6 x 35.4
lass (without oil)	kg Ib	80 176	140 309	160 353
Power supply	-	230 V ±10%, 50/60 Hz Max 3000 W	230 V ±10%, 50/60 Hz Max 3000 W	230 V ±10%, 50/60 Hz Max 3000 W
Data connections	-	Ethernet (for PC) USB (for printer)	Ethernet (for PC) USB (for printer)	Ethernet (for PC) USB (for printer)

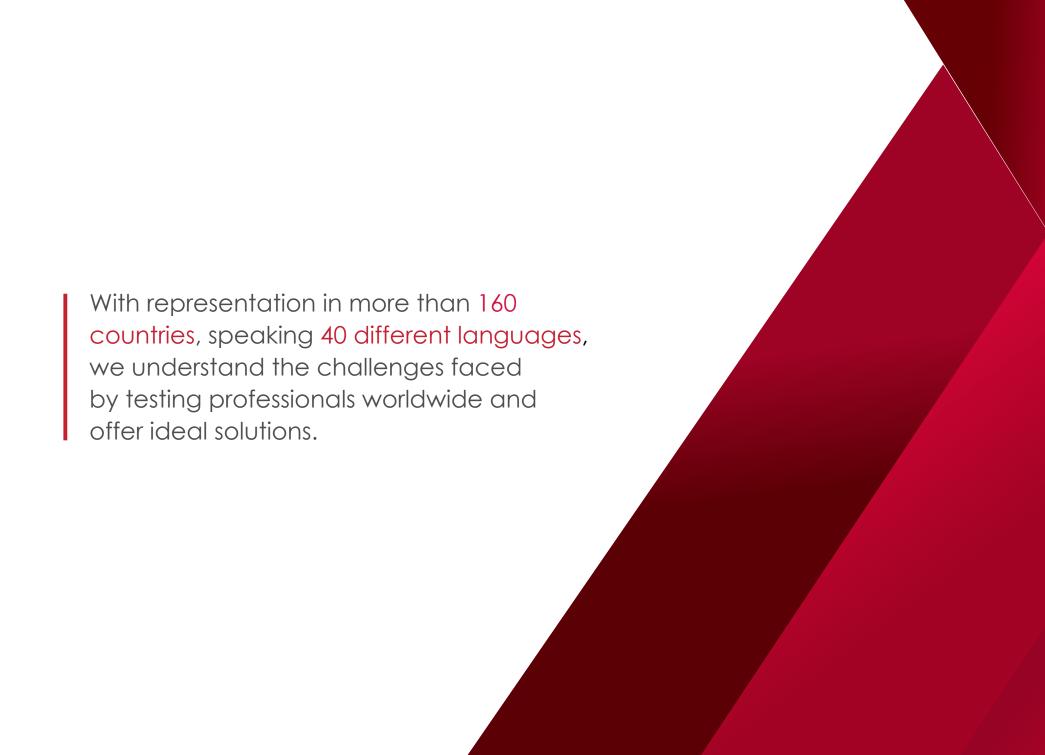



"True innovation occurs when our product designers and developers show relentless curiosity towards the needs of our customers. This builds an understanding that allows us to anticipate and create a new suite of solutions that are Simpler, Smarter, and Safer."

Yahya Gharagozlou

Group President ITW Test & Measurement

www.instron.com

Instron is a registered trademark of Illinois Tool Works Inc. (ITW). Other names, logos, icons and marks identifying Instron products and services referenced herein are trademarks of ITW and may not be used without the prior written permission of ITW. Other product and company names listed are trademarks or trade names of their respective companies. Copyright © 2016 Illinois Tool Works Inc. All rights reserved. All of the specifications shown in this document are subject to change without notice.